

INSTRUCTIONS FOR USE

- IMPORTANT!** Read and follow all instructions before use (contact reseller for details or to obtain a copy).
- TLE (Top Lip Extension) is an accessory for some EZ-ACCESS ramps. **NOTE:** When using with SUITCASE ramp, refer to WARNINGS and steps 1-7 only.
- Ensure that you have received two each of the TLE (they must be installed in pairs).
- Open ramp (refer to ramp instructions) and insert one TLE under the TOP transition plate as shown in FIGURE 1. Rotate the TLE upward as shown in FIGURE 2.
- Ensure that the TLE is firmly seated in the TOP transition plate and that the square hole in the TOP transition plate is aligned with the hole in the TLE.
- Install the carriage bolt from the topside of TLE, through the square hole. Secure with washer and nut on underside of TLE as shown in FIGURE 3 and tighten the nut using a 7/16" wrench.
- Repeat for other side of ramp.
- When using with TRIFOLD ramp, install the hook and loop strap onto one TLE as shown in FIGURE 4.
- Fold ramp to storage position (refer to ramp instructions), then secure strap to the other TLE for transport, etc. as shown in FIGURE 5.

- ▲ **WARNING!** DO NOT USE ramp with TLE unless the carriage bolt and nut are secure (see FIGURE 3).
- ▲ **WARNING!** Ensure that the top plate is secure before using. We recommend attaching the ramp temporarily to the landing using the provided clevis pins (drill 5/16" hole in landing) or semi-permanently using screws through the pre-drilled holes at the top of each section. To prevent ramp from becoming dislodged during use, make sure that the top plate overlaps the resting surface as far as possible. The surface holding the top plate needs to be substantial; make sure that the surface and supporting structure is of sound construction.
- ▲ **WARNING!** Ensure the TLE's surface area is resting fully on the landing during use (see FIGURE 6).
- ▲ **WARNING!** TLE is designed for use with unoccupied mobility equipment only.
- ▲ **WARNING!** Use of the TLE reduces overall ramp weight capacity to 500 lbs.

SAFETY & MAINTENANCE

Regularly check ramp for damage, and check and tighten any loose fasteners. Keep ramp surface clear of ice, snow, dirt, leaves, and other debris. For snow/ice, remove accumulation and sweep the surface clean. To remove dirt, leaves, and other debris, sweeping the surface or using a garden hose will usually suffice, but, if needed, use a damp cloth or soft brush with non-alkaloid detergent and water. Rinse well and use extra caution if ramp surface is wet.

FIGURE 1

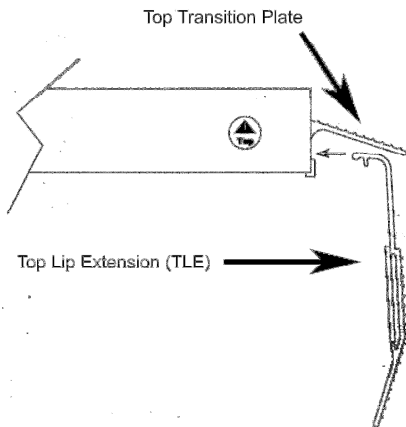


FIGURE 2

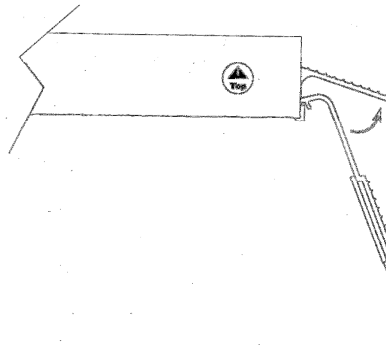


FIGURE 3

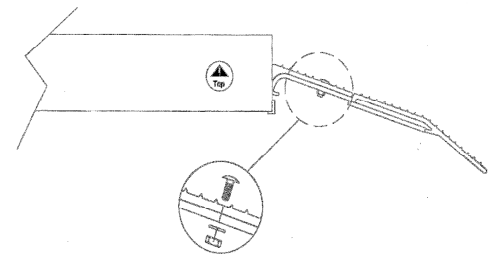


FIGURE 4

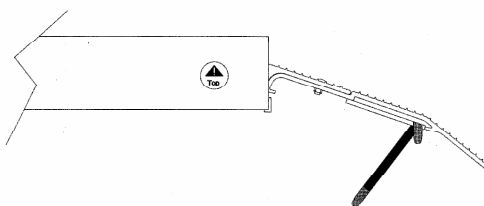


FIGURE 5

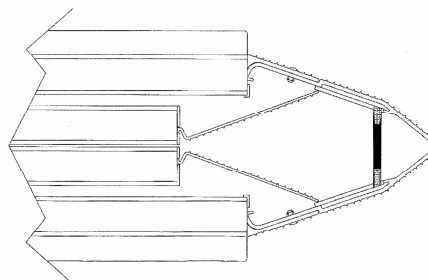


FIGURE 6

